
Email : sales@scilabexport.com
Phone: +91-7082934803

Product Name :

Double Gel Electrophoresis Apparatus

Product Code :

SLE/BL/29

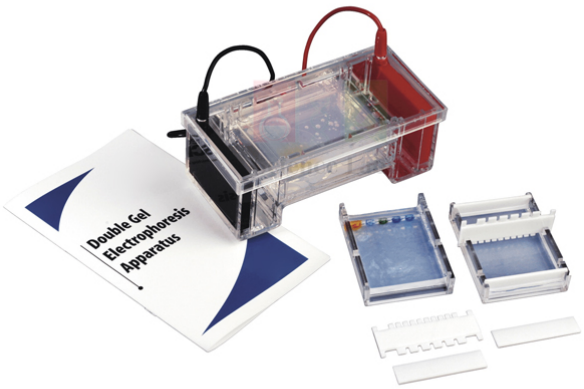

Science Lab

Email : sales@scilabexport.com
Phone: +91-7082934803

Description :

The **Double Gel Electrophoresis Apparatus** (Product Code: SLE/BL/29) is a horizontal gel tank that runs two gels at the same time, so two student groups can share one unit or a class can compare two sets of samples in a single session. Its leak-free, one-piece moulded chamber is made of durable clear acrylic and fitted with pure platinum electrodes, and a lid interlock prevents use until the lid is locked in place. Schools, colleges and biotechnology teaching labs use it as a complete gel electrophoresis kit for DNA and dye separation practicals. It is made by a biologylab equipment manufacturer in India, **Jain Scientific Equipments Pvt Ltd (ScienceLab)**,Ambala.

Specifications

Parameter	Detail
Capacity	Runs two gels simultaneously
Chamber	Leak-free, one-piece moulded, durable clear acrylic
Electrodes	Pure platinum
Casting trays	Two; hold melted agarose without tape or snap-on ends
Combs	Six 12-well forming combs
Well position	At the end or the centre of the gel
Connections	Patch cords included
Safety	Chamber cannot be used unless the lid is locked in place
Tank size	155 mm (L) x 88 mm (W) x 62.5 mm (H)
Tray size	100 mm (L) x 75 mm (W)
Power supply	Not included; available separately

Typical Use

Students place a comb in one of the specially designed casting trays and pour melted agarose straight in – no taping of the tray ends is needed. Setting the comb at the end of the gel suits DNA, which moves in one direction, while a centre position lets charged dyes of opposite sign migrate both ways from the wells. Once set, the combs are lifted out, the trays go into the tank, running buffer is added, and samples are loaded into the 12 wells of each gel with a micropipette.

With the lid locked and the patch cords connected to a suitable power supply, both gels run together, which halves the time needed to get a whole class through the practical. In a biotechnology lab kit setting this makes it practical to compare, for example, a sample and a reference side by side. Rinse the chamber with water after use, avoid solvents and abrasive cleaners on the acrylic, and handle the platinum electrode wires gently when cleaning.

Manufacturer

Scien ana
manu
equip
more