Email : sales@scilabexport.com
Phone: +91-7082934803


**Product Name :**

Prism – Acrylic, Set of 7

**Product Code :**

SLE/PHY/274

Image not found or type unknown

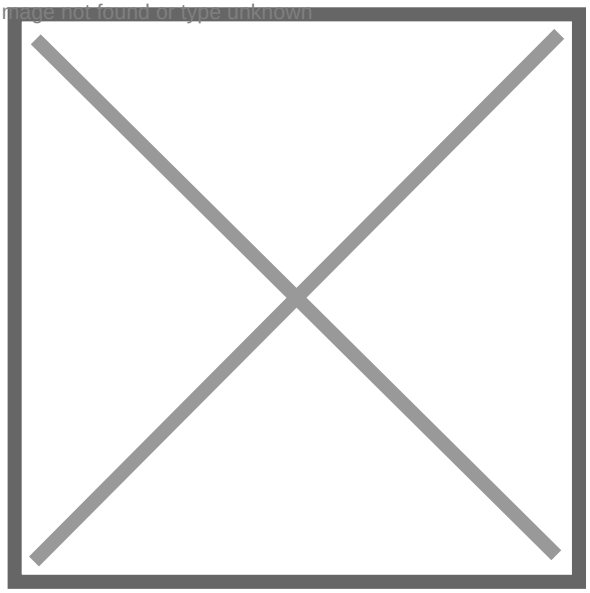

Email : sales@scilabexport.com
Phone: +91-7082934803

## Description :

**Prism – Acrylic, Set of 7** (Product Code: SLE/PHY/274) is a set of seven polished acrylic optical blocks used to study reflection, refraction and total internal reflection with a ray box in school and college optics lessons. Every block is 25 mm thick with all faces fully polished, and the shapes cover bi-convex, bi-concave, rectangular, semi-circular and three triangular prisms. It is made by **Jain Scientific Equipments Pvt Ltd (ScienceLab)**, an OEM manufacturer, at its own factory in Ambala, India.

### Specifications

| Parameter | Detail |
|---|---|
| Number of blocks | 7 |
| Material | Acrylic |
| Thickness | 25 mm (all blocks) |
| Finish | All faces fully polished |
| Bi-convex block | 75 mm long curved faces, 100 mm radius |
| Bi-concave block | 75 mm long curved faces, 100 mm radius |
| Rectangular block | 75 x 50 mm |
| Semi-circular block | 75 mm diameter |
| Equilateral triangle | 60° x 60° x 60°, 58 mm side |
| Right-angled isosceles triangle | 90° x 45° x 45°, 75 mm hypotenuse |
| 30-60-90 triangle | 90° x 60° x 30°, 75 mm hypotenuse |
| Product code | SLE/PHY/274 |
| Manufacturer | Jain Scientific Equipments Pvt Ltd (ScienceLab) — OEM |
| Country of origin | India (made in Ambala, Haryana) |

### What Is the Acrylic Prism Set Used For?

Laid flat on paper in the path of a ray box beam, each block demonstrates a different part of the optics syllabus. The rectangular block shows lateral displacement and lets students work out refractive index from the angles of incidence and refraction; the semi-circular block is the standard shape for finding the critical angle and showing total internal reflection; and the bi-convex and bi-concave blocks make parallel rays converge or diverge like thin lenses drawn in two dimensions.

The 60° triangle disperses white light into a spectrum, while the two right-angled triangles show how a prism can turn a ray through a right angle or send it back the way it came by total internal reflection, the principle behind periscopes and binoculars. As a compact optics experiment kit for schoo r, and c

Wha